

TOPO TA Cloning[®]

Version R
8 April 2004
25-0184

TOPO TA Cloning[®]

Five-minute cloning of *Taq* polymerase-amplified PCR products

Catalog nos. K4500-01, K4500-40, K4510-20, K4520-01, K4520-40, K4550-01,
K4550-40, K4560-01, K4560-40 (pCR[®] 2.1-TOPO[®])

Catalog nos. K4600-01, K4600-40, K4610-20, K4620-01, K4620-40, K4650-01,
K4650-40, K4660-01, K4660-40 (pCR[®] II-TOPO[®])

A Limited Label License covers this product (see Purchaser Notification). By use of this product, you accept the terms and conditions of the Limited Label License.



www.invitrogen.com
tech_service@invitrogen.com

Table of Contents

Table of Contents	iii
Kit Contents and Storage	v
Accessory Products.....	viii
Methods	1
Overview.....	1
Producing PCR Products	3
Setting Up the TOPO [®] Cloning Reaction.....	4
General Guidelines for Transforming Competent Cells	6
Transforming One Shot [®] Mach1 [™] -T1 ^R Competent Cells	7
Transforming One Shot [®] DH5 α [™] -T1 ^R , TOP10, and TOP10F [′] Competent Cells.....	9
Analyzing Transformants.....	12
Optimizing the TOPO [®] Cloning Reaction	14
Map of pCR [®] 2.1-TOPO [®]	15
Map of pCR [®] II-TOPO [®]	16
Performing the Control Reactions	17
Appendix	20
Purifying PCR Products.....	20
Addition of 3′ A-Overhangs Post-Amplification	22
Recipes.....	23
Technical Service.....	24
Product Qualification	26
Purchaser Notification	27
References	30

Kit Contents and Storage

Shipping and Storage

TOPO TA Cloning[®] Kits are shipped on dry ice. Each kit contains a box with TOPO TA Cloning[®] reagents (Box 1) and a box with One Shot[®] Chemically Competent or Electrocomp[™] cells (Box 2).

Store Box 1 at -20°C and Box 2 at -80°C.

Types of TOPO TA Cloning[®] Kits

TOPO TA Cloning[®] Kits are available with either pCR[®]2.1-TOPO[®] or pCR[®]II-TOPO[®] and either DH5 α [™]-T1^R, Mach1[™]-T1^R, TOP10, or TOP10F[′] One Shot[®] Chemically Competent cells or TOP10 One Shot[®] Electrocomp[™] cells (see page vii for the genotypes of the strains).

Product	Reactions	One Shot [®] Cells	Type of Cells	Catalog no.
TOPO TA Cloning [®] Kit (containing pCR [®] 2.1-TOPO [®])	20	TOP10	chem. competent	K4500-01
	40	TOP10	chem. competent	K4500-40
	20	Mach1 [™] -T1 ^R	chem. competent	K4510-20
	20	DH5 α [™] -T1 ^R	chem. competent	K4520-01
	40	DH5 α [™] -T1 ^R	chem. competent	K4520-40
	20	TOP10F [′]	chem. competent	K4550-01
	40	TOP10F [′]	chem. competent	K4550-40
	20	TOP10	electrocompetent	K4560-01
	40	TOP10	electrocompetent	K4560-40
TOPO TA Cloning [®] Kit Dual Promoter (containing pCR [®] II-TOPO [®])	20	TOP10	chem. competent	K4600-01
	40	TOP10	chem. competent	K4600-40
	20	Mach1 [™] -T1 ^R	chem. competent	K4610-20
	20	DH5 α [™] -T1 ^R	chem. competent	K4620-01
	40	DH5 α [™] -T1 ^R	chem. competent	K4620-40
	20	TOP10F [′]	chem. competent	K4650-01
	40	TOP10F [′]	chem. competent	K4650-40
	20	TOP10	electrocompetent	K4660-01
	40	TOP10	electrocompetent	K4660-40

continued on next page

Kit Contents and Storage, continued

TOPO TA Cloning® Reagents

TOPO TA Cloning® reagents (Box 1) are listed below. Note that the user must supply Taq polymerase. Store Box 1 at -20°C.

Item	Concentration	Amount
pCR®2.1-TOPO® or pCR®II-TOPO®	10 ng/μl plasmid DNA in: 50% glycerol 50 mM Tris-HCl, pH 7.4 (at 25°C) 1 mM EDTA 1 mM DTT 0.1% Triton X-100 100 μg/ml BSA phenol red	20 μl
10X PCR Buffer	100 mM Tris-HCl, pH 8.3 (at 42°C) 500 mM KCl 25 mM MgCl ₂ 0.01% gelatin	100 μl
Salt Solution	1.2 M NaCl 0.06 M MgCl ₂	50 μl
dNTP Mix	12.5 mM dATP; 12.5 mM dCTP 12.5 mM dGTP; 12.5 mM dTTP neutralized at pH 8.0 in water	10 μl
M13 Forward (-20) Primer	0.1 μg/μl in TE Buffer	20 μl
M13 Reverse Primer	0.1 μg/μl in TE Buffer	20 μl
Control Template	0.1 μg/μl in TE Buffer	10 μl
Control PCR Primers	0.1 μg/μl each in TE Buffer	10 μl
Sterile Water	--	1 ml

Sequence of Primers

The table below describes the sequence and pmoles supplied of the sequencing primers included in this kit.

Primer	Sequence	pMoles Supplied
M13 Forward (-20)	5'-GTAAAACGACGGCCAG-3'	407
M13 Reverse	5'-CAGGAAACAGCTATGAC-3'	385

continued on next page

Kit Contents and Storage, continued

One Shot® Reagents

The table below describes the items included in each One Shot® competent cell kit. Store at -80°C.

Item	Composition	Amount
S.O.C. Medium (may be stored at +4°C or room temperature)	2% Tryptone 0.5% Yeast Extract 10 mM NaCl 2.5 mM KCl 10 mM MgCl ₂ 10 mM MgSO ₄ 20 mM glucose	6 ml
TOP10, Mach1™-T1 ^R , DH5α™-T1 ^R , or TOP10F' cells OR TOP10 cells	Chemically Competent Electrocomp™	21 x 50 µl
pUC19 Control DNA	10 pg/µl in 5 mM Tris-HCl, 0.5 mM EDTA, pH 8	50 µl

Genotypes of *E. coli* Strains

DH5α™-T1^R: Use this strain for general cloning and blue/white screening without IPTG. Strain is resistant to T1 bacteriophage.

F⁻ ϕ 80*lacZ*ΔM15 Δ(*lacZYA-argF*)U169 *recA1 endA1 hsdR17*(r_k⁻, m_k⁺) *phoA supE44 thi-1 gyrA96 relA1 tonA* (confers resistance to phage T1)

Mach1™-T1^R: Use this strain for general cloning and blue/white screening without IPTG. Strain is resistant to T1 bacteriophage.

F⁻ ϕ 80(*lacZ*)ΔM15 Δ*lacX74 hsdR*(r_k⁻, m_k⁺) Δ*recA1398 endA1 tonA* (confers resistance to phage T1)

TOP10: Use this strain for general cloning and blue/white screening without IPTG.

F⁻ *mcrA* Δ(*mrr-hsdRMS-mcrBC*) Φ 80*lacZ*ΔM15 Δ*lacX74 recA1 araD139* Δ(*ara-leu*)7697 *galU galK rpsL* (Str^R) *endA1 nupG*

TOP10F': This strain overexpresses the Lac repressor (*lacI^q* gene). For blue/white screening, you will need to add IPTG to the plates to obtain expression from the *lac* promoter. This strain contains the F episome and can be used for single-strand rescue of plasmid DNA containing an f1 origin.

F' {*lacI^q Tn10* (Tet^R)} *mcrA* Δ(*mrr-hsdRMS-mcrBC*) Φ 80*lacZ*ΔM15 Δ*lacX74 recA1 araD139* Δ(*ara-leu*)7697 *galU galK rpsL* (Str^R) *endA1 nupG*

Information for Non-U.S. Customers Using Mach1™-T1^R Cells

The parental strain of Mach1™-T1^R *E. coli* is the non-K-12, wild-type W strain (ATCC #9637, S. A. Waksman). Although the parental strain is generally classified as Biosafety Level 1 (BL-1), we recommend that you consult the safety department of your institution to verify the Biosafety Level.

Accessory Products

Additional Products

The table below lists additional products that may be used with TOPO® TA Cloning Kits. For more information, refer to our Web site (www.invitrogen.com) or contact Technical Service (page 24).

Item	Amount	Catalog no.
<i>Taq</i> DNA Polymerase, Native	100 units	18038-018
	500 units	18038-042
<i>Taq</i> DNA Polymerase, Recombinant	100 units	10342-053
	500 units	10342-020
Platinum® <i>Taq</i> DNA Polymerase High Fidelity	100 units	11304-011
One Shot® TOP10 Chemically Competent <i>E. coli</i>	10 reactions	C4040-10
	20 reactions	C4040-03
	40 reactions	C4040-06
One Shot® TOP10 Electrocompetent <i>E. coli</i>	10 reactions	C4040-50
	20 reactions	C4040-52
One Shot® Mach1™-T1 ^R Chemically Competent <i>E. coli</i>	20 reactions	C8620-03
One Shot® MAX Efficiency® DH5α-T1 ^R Chemically Competent <i>E. coli</i>	20 reactions	12297-016
One Shot® TOP10F' Chemically Competent <i>E. coli</i>	20 reactions	C3030-03
	40 reactions	C3030-06
S.N.A.P.™ MidiPrep Kit	20 reactions	K1910-01
Ampicillin	200 mg	11593-019
Kanamycin	5 g	11815-024
	25 g	11815-032
	100 ml (10 mg/ml)	18160-054
X-gal	100 mg	15520-034
	1 g	15520-018
IPTG	1 g	15529-019
S.O.C. Medium	10 x 10 ml	15544-034

Methods

Overview

Introduction

TOPO TA Cloning[®] provides a highly efficient, 5-minute, one-step cloning strategy ("TOPO[®] Cloning") for the direct insertion of *Taq* polymerase-amplified PCR products into a plasmid vector. No ligase, post-PCR procedures, or PCR primers containing specific sequences are required.


How It Works

The plasmid vector (pCR[®]II-TOPO[®] or pCR[®]2.1-TOPO[®]) is supplied linearized with:

- Single 3'-thymidine (T) overhangs for TA Cloning[®]
- Topoisomerase I covalently bound to the vector (referred to as "activated" vector)

Taq polymerase has a nontemplate-dependent terminal transferase activity that adds a single deoxyadenosine (A) to the 3' ends of PCR products. The linearized vector supplied in this kit has single, overhanging 3' deoxythymidine (T) residues. This allows PCR inserts to ligate efficiently with the vector.

Topoisomerase I from *Vaccinia* virus binds to duplex DNA at specific sites and cleaves the phosphodiester backbone after 5'-CCCTT in one strand (Shuman, 1991). The energy from the broken phosphodiester backbone is conserved by formation of a covalent bond between the 3' phosphate of the cleaved strand and a tyrosyl residue (Tyr-274) of topoisomerase I. The phospho-tyrosyl bond between the DNA and enzyme can subsequently be attacked by the 5' hydroxyl of the original cleaved strand, reversing the reaction and releasing topoisomerase (Shuman, 1994).

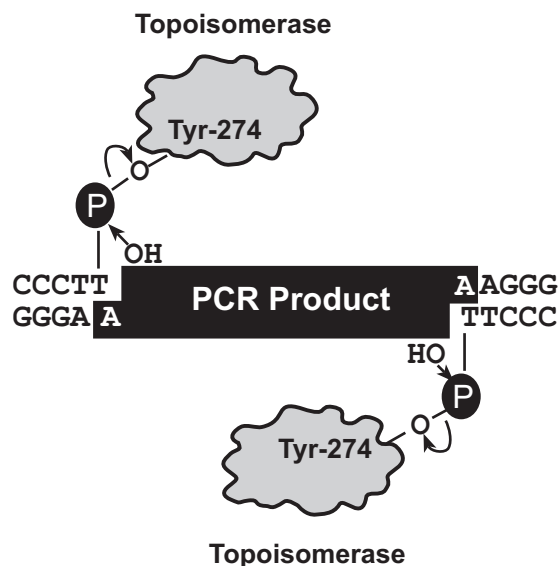

continued on next page

Overview, continued

Experimental Outline

- Produce Your PCR Product
 - Set Up TOPO[®] Cloning Reaction (Mix Together PCR Product and TOPO[®] Vector)
 - Incubate 5 Minutes at Room Temperature
 - Transform TOPO[®] Cloning Reaction into One Shot[®] Competent Cells
 - Select and Analyze 10 White or Light Blue Colonies for Insert
-

Producing PCR Products

Introduction

It is important to properly design your PCR primers to ensure that you obtain the product you need for your studies. Once you have decided on a PCR strategy and have synthesized the primers, you are ready to produce your PCR product. **Remember that your PCR product will have single 3' adenine overhangs.**



Note

Do not add 5' phosphates to your primers for PCR. The PCR product synthesized will not ligate into pCR[®]2.1-TOPO[®] or pCR[®]II-TOPO[®].

Materials Supplied by the User

You will need the following reagents and equipment.

- *Taq* polymerase
 - Thermocycler
 - DNA template and primers for PCR product
-

Polymerase Mixtures

If you wish to use a mixture containing *Taq* polymerase and a proofreading polymerase, *Taq* must be used in excess of a 10:1 ratio to ensure the presence of 3' A-overhangs on the PCR product.

If you use polymerase mixtures that do not have enough *Taq* polymerase or a proofreading polymerase only, you can add 3' A-overhangs using the method on page 22.

Producing PCR Products

1. Set up the following 50 μ l PCR reaction. Use less DNA if you are using plasmid DNA as a template and more DNA if you are using genomic DNA as a template. Use the cycling parameters suitable for your primers and template. Be sure to include a 7 to 30 minute extension at 72°C after the last cycle to ensure that all PCR products are full length and 3' adenylated.

DNA Template	10-100 ng
10X PCR Buffer	5 μ l
50 mM dNTPs	0.5 μ l
Primers (100-200 ng each)	1 μ M each
Sterile water	add to a final volume of 49 μ l
<i>Taq</i> Polymerase (1 unit/ μ l)	1 μ l
Total Volume	50 μ l

2. Check the PCR product by agarose gel electrophoresis. You should see a single, discrete band. If you do not see a single band, refer to the **Note** below.
-



Note

If you do not obtain a single, discrete band from your PCR, you may gel-purify your fragment before using the TOPO TA Cloning[®] Kit (see page 20). Take special care to avoid sources of nuclease contamination. Alternatively, you may optimize your PCR to eliminate multiple bands and smearing (Innis *et al.*, 1990). The PCR Optimizer[™] Kit (Catalog no. K1220-01) incorporates many of the recommendations found in this reference. For more information, refer to our Web site (www.invitrogen.com) or contact Technical Service (page 24).

Setting Up the TOPO[®] Cloning Reaction

Introduction

Once you have produced the desired PCR product, you are ready to TOPO[®] Clone it into the pCR[®]2.1-TOPO[®] or pCR[®]II-TOPO[®] vector and transform the recombinant vector into competent *E. coli*. It is important to have everything you need set up and ready to use to ensure that you obtain the best possible results. We suggest that you read this section and the sections detailing transformation of competent cells (pages 6-10) before beginning. If this is the first time you have TOPO[®] Cloned, perform the control reactions on pages 17-18 in parallel with your samples.



Note

Recent experiments at Invitrogen demonstrate that inclusion of salt (200 mM NaCl; 10 mM MgCl₂) in the TOPO[®] Cloning reaction increases the number of transformants 2- to 3-fold. We have also observed that in the presence of salt, incubation times of greater than 5 minutes can also increase the number of transformants. This is in contrast to earlier experiments **without salt** where the number of transformants decreases as the incubation time increases beyond 5 minutes.

Inclusion of salt allows for longer incubation times because it prevents topoisomerase I from rebinding and potentially nicking the DNA after ligating the PCR product and dissociating from the DNA. The result is more intact molecules leading to higher transformation efficiencies.



Important

Because of the above results, we recommend adding salt to the TOPO[®] Cloning reaction. A stock salt solution is provided in the kit for this purpose. **Note that the amount of salt added to the TOPO[®] Cloning reaction varies depending on whether you plan to transform chemically competent cells or electrocompetent cells (see below).** For this reason two different TOPO[®] Cloning reactions are provided to help you obtain the best possible results. Read the following information carefully.

Transforming Chemically Competent *E. coli*

For TOPO[®] Cloning and transformation into chemically competent *E. coli*, adding sodium chloride and magnesium chloride to a final concentration of 200 mM NaCl, 10 mM MgCl₂ in the TOPO[®] Cloning reaction increases the number of colonies over time. A Salt Solution (1.2 M NaCl; 0.06 M MgCl₂) is provided to adjust the TOPO[®] Cloning reaction to the recommended concentration of NaCl and MgCl₂.

Transforming Electrocompetent *E. coli*

For TOPO[®] Cloning and transformation of electrocompetent *E. coli*, salt must also be included in the TOPO[®] Cloning reaction, but the amount of salt **must be reduced** to 50 mM NaCl, 2.5 mM MgCl₂ to prevent arcing. The Salt Solution is diluted 4-fold to prepare a 300 mM NaCl, 15 mM MgCl₂ solution for convenient addition to the TOPO[®] Cloning reaction (see next page).

continued on next page

Setting Up the TOPO[®] Cloning Reaction, continued

Setting Up the TOPO[®] Cloning Reaction

The table below describes how to set up your TOPO[®] Cloning reaction (6 μ l) for eventual transformation into either chemically competent or electrocompetent TOP10 or chemically competent DH5 α [™]-T1^R, Mach1[™]-T1^R, or TOP10F' One Shot[®] *E. coli*. Additional information on optimizing the TOPO[®] Cloning reaction for your needs can be found on page 14.

Note: The red color of the TOPO[®] vector solution is normal and is used to visualize the solution.

Reagent*	Chemically Competent <i>E. coli</i>	Electrocompetent <i>E. coli</i>
Fresh PCR product	0.5 to 4 μ l	0.5 to 4 μ l
Salt Solution	1 μ l	--
Dilute Salt Solution	--	1 μ l
Sterile Water	add to a total volume of 5 μ l	add to a total volume of 5 μ l
TOPO [®] vector	1 μ l	1 μ l
Final Volume	6 μl	6 μl

*Store all reagents at -20°C when finished. Salt solutions and water can be stored at room temperature or +4°C.

Performing the TOPO[®] Cloning Reaction

1. Mix reaction gently and incubate for **5 minutes** at room temperature (22-23°C).
Note: For most applications, 5 minutes will yield plenty of colonies for analysis. Depending on your needs, the length of the TOPO[®] Cloning reaction can be varied from 30 seconds to 30 minutes. For routine subcloning of PCR products, 30 seconds may be sufficient. For large PCR products (> 1 kb) or if you are TOPO[®] Cloning a pool of PCR products, increasing the reaction time will yield more colonies.
2. Place the reaction on ice and proceed to **General Guidelines for Transforming Competent Cells**, next page.

Note: You may store the TOPO[®] Cloning reaction at -20°C overnight.

General Guidelines for Transforming Competent Cells

Introduction

Once you have performed the TOPO[®] Cloning reaction, you will transform your pCR[®]2.1-TOPO[®] or pCR[®]II-TOPO[®] construct into the competent *E. coli* provided with your kit. General guidelines for transformation are provided below. For transformation protocols, refer to the section entitled **Transforming One Shot[®] Mach1[™]-T1^R Competent Cells** (pages 7-8) or **Transforming One Shot[®] DH5 α [™]-T1^R, TOP10, and TOP10F[′] Competent Cells** (pages 9-11) depending on the competent *E. coli* you wish to transform.

Selecting a One Shot[®] Chemical Transformation Protocol

Two protocols are provided to transform One Shot[®] chemically competent *E. coli*. Consider the following factors when choosing the protocol that best suits your needs.

If you wish to...	Then use the...
maximize the number of transformants	regular chemical transformation protocol
clone large PCR products (>1000 bp)	
use kanamycin as the selective agent (see Important note below)	
obtain transformants as quickly as possible	rapid chemical transformation protocol



Important

If you will be using kanamycin as the selective agent for chemical transformation, use the regular chemical transformation protocol. The rapid chemical transformation protocol is only suitable for transformations using ampicillin selection.



If you use a plasmid template for your PCR that carries either the ampicillin or kanamycin resistance marker, we recommend that you use the other selection agent to select for transformants. For example, if the plasmid template contains the ampicillin resistance marker, then use kanamycin to select for transformants. The template is carried over into the TOPO[®] Cloning and transformation reactions, resulting in transformants that are ampicillin resistant and white, but are not the desired construct.

Transforming One Shot[®] Mach1[™]-T1^R Competent Cells

Introduction

Protocols to transform One Shot[®] Mach1[™]-T1^R chemically competent *E. coli* are provided below. If are transforming cells other than Mach1[™]-T1^R cells, refer to the section entitled **Transforming One Shot[®] DH5 α [™]-T1^R, TOP10, and TOP10F['] Competent Cells** (pages 9-11).



Note

The Mach1[™]-T1^R strain allows you to visualize colonies 8 hours after plating on ampicillin selective plates. If you are using kanamycin selection, you will need to incubate plates overnight in order to visualize colonies.

With the Mach1[™]-T1^R strain, you may also prepare plasmid DNA 4 hours after inoculating a single, overnight-grown colony. Note that you will get sufficient growth of transformed cells within 4 hours in either ampicillin or kanamycin selective media.

Materials Supplied by the User

In addition to general microbiological supplies (*e.g.* plates, spreaders), you will need the following reagents and equipment.

- TOPO[®] Cloning reaction from Performing the TOPO[®] Cloning Reaction, Step 2 (page 5)
 - S.O.C. medium (included with the kit)
 - LB plates containing 50 μ g/ml ampicillin or 50 μ g/ml kanamycin
 - 40 mg/ml X-gal in dimethylformamide (DMF)
 - 42°C water bath
 - 37°C shaking and non-shaking incubator
-

Preparing for Transformation

For each transformation, you will need one vial of competent cells and two selective plates.

- Equilibrate a water bath to 42°C.
 - Warm the vial of S.O.C. medium from Box 2 to room temperature.
 - Warm selective plates at 37°C for 30 minutes (see Important note below).
 - Spread 40 μ l of 40 mg/ml X-gal on each LB plate and incubate at 37°C until ready for use.
 - Thaw **on ice** 1 vial of One Shot[®] cells for each transformation.
-



Important

If you are performing the rapid chemical transformation protocol or if you wish to visualize colonies within 8 hours of plating, it is essential that you prewarm your LB plates containing 50-100 μ g/ml ampicillin prior to spreading.

continued on next page

Transforming One Shot[®] Mach1[™]-T1^R Competent Cells, continued

One Shot[®] Chemical Transformation Protocol

For optimal growth of Mach1[™]-T1^R *E. coli* cells, it is essential that selective plates are prewarmed to 37°C prior to spreading.

1. Add 2 µl of the TOPO[®] Cloning reaction from Performing the TOPO[®] Cloning Reaction, Step 2, page 5 into a vial of One Shot[®] Chemically Competent *E. coli* and mix gently. **Do not mix by pipetting up and down.**
2. Incubate on ice for 5 to 30 minutes.
Note: Longer incubations on ice do not seem to have any affect on transformation efficiency. The length of the incubation is at the user's discretion.
3. Heat-shock the cells for 30 seconds at 42°C without shaking.
4. Immediately transfer the tubes to ice.
5. Add 250 µl of room temperature S.O.C. medium.
6. Cap the tube tightly and shake the tube horizontally (200 rpm) at 37°C for 1 hour.
7. Spread 10-50 µl from each transformation on a **prewarmed** selective plate. To ensure even spreading of small volumes, add 20 µl of S.O.C. medium. We recommend that you plate two different volumes to ensure that at least one plate will have well-spaced colonies.
8. Incubate plates at 37°C. If you are using ampicillin selection, visible colonies should appear within 8 hours, and blue/white screening can be performed after 12 hours. For kanamycin selection, incubate plates overnight.
9. An efficient TOPO[®] Cloning reaction should produce several hundred colonies. Pick ~10 white or light blue colonies for analysis (see Analyzing Positive Clones, page 12). Do not pick dark blue colonies.

Rapid One Shot[®] Chemical Transformation Protocol

An alternative protocol is provided below for rapid transformation of One Shot[®] Mach1[™]-T1^R cells. This protocol is **only** recommended for transformations using **ampicillin** selection. For more information on selecting a transformation protocol, refer to page 6.

Note: It is essential that LB plates containing ampicillin are prewarmed to 37°C prior to spreading.

1. Add 4 µl of the TOPO[®] Cloning reaction from Performing the TOPO[®] Cloning Reaction, Step 2, page 5 into a vial of One Shot[®] Chemically Competent *E. coli* and mix gently. **Do not mix by pipetting up and down.**
 2. Incubate on ice for 5 minutes.
 3. Spread 50 µl of cells on a prewarmed LB plate containing 50-100 µg/ml ampicillin and incubate overnight at 37°C.
 4. An efficient TOPO[®] Cloning reaction should produce several hundred colonies. Pick ~10 white or light blue colonies for analysis (see Analyzing Positive Clones, page 12). Do not pick dark blue colonies.
-

Transforming One Shot[®] DH5 α [™]-T1^R, TOP10, and TOP10F['] Competent Cells

Introduction

Protocols to transform One Shot[®] DH5 α [™]-T1^R, TOP10, and TOP10F['] competent *E. coli* are provided below. Both chemical transformation and electroporation protocols are provided. If you are transforming Mach1[™]-T1^R cells, refer to the section entitled **Transforming One Shot[®] Mach1[™]-T1^R Competent Cells** (pages 7-8).

Materials Supplied by the User

In addition to general microbiological supplies (*e.g.* plates, spreaders), you will need the following reagents and equipment.

- TOPO[®] Cloning reaction from Performing the TOPO[®] Cloning Reaction, Step 2 (page 5)
 - S.O.C. medium (included with the kit)
 - LB plates containing 50 μ g/ml ampicillin or 50 μ g/ml kanamycin
 - 40 mg/ml X-gal in dimethylformamide (DMF)
 - 100 mM IPTG in water (for use with TOP10F['])
 - 15 ml snap-cap plastic culture tubes (sterile) (electroporation only)
 - 42°C water bath or an electroporator and 0.1 or 0.2 cm cuvettes
 - 37°C shaking and non-shaking incubator
-

Preparation for Transformation

For each transformation, you will need one vial of competent cells and two selective plates.

- Equilibrate a water bath to 42°C (for chemical transformation) or set up your electroporator.
 - Warm the vial of S.O.C. medium from Box 2 to room temperature.
 - Warm selective plates at 37°C for 30 minutes (see Important note below).
 - Spread 40 μ l of 40 mg/ml X-gal on each LB plate and incubate at 37°C until ready for use.
 - For TOP10F['] cells, spread 40 μ l of 100 mM IPTG in addition to X-gal on each LB plate and incubate at 37°C until ready for use. IPTG is required for blue/white screening.
 - Thaw **on ice** 1 vial of One Shot[®] cells for each transformation.
-



Important

If you are performing the rapid chemical transformation protocol, it is essential that you prewarm your LB plates containing 50-100 μ g/ml ampicillin prior to spreading.

continued on next page

Transforming One Shot[®] DH5 α [™] -T1^R, TOP10, and TOP10F['] Competent Cells, continued

One Shot[®] Chemical Transformation Protocol

1. Add 2 μ l of the TOPO[®] Cloning reaction from Performing the TOPO[®] Cloning Reaction, Step 2, page 5 into a vial of One Shot[®] Chemically Competent *E. coli* and mix gently. **Do not mix by pipetting up and down.**
 2. Incubate on ice for 5 to 30 minutes.
Note: Longer incubations on ice do not seem to have any affect on transformation efficiency. The length of the incubation is at the user's discretion.
 3. Heat-shock the cells for 30 seconds at 42°C without shaking.
 4. Immediately transfer the tubes to ice.
 5. Add 250 μ l of room temperature S.O.C. medium.
 6. Cap the tube tightly and shake the tube horizontally (200 rpm) at 37°C for 1 hour.
 7. Spread 10-50 μ l from each transformation on a prewarmed selective plate and incubate overnight at 37°C. To ensure even spreading of small volumes, add 20 μ l of S.O.C. medium. We recommend that you plate two different volumes to ensure that at least one plate will have well-spaced colonies.
 8. An efficient TOPO[®] Cloning reaction should produce several hundred colonies. Pick ~10 white or light blue colonies for analysis (see Analyzing Positive Clones, page 12). Do not pick dark blue colonies.
-

Rapid One Shot[®] Chemical Transformation Protocol

An alternative protocol is provided below for rapid transformation of One Shot[®] chemically competent *E. coli*. This protocol is **only** recommended for transformations using **ampicillin** selection. For more information on selecting a transformation protocol, refer to page 6.

Note: It is essential that LB plates containing ampicillin are prewarmed prior to spreading.

1. Add 4 μ l of the TOPO[®] Cloning reaction from Performing the TOPO[®] Cloning Reaction, Step 2, page 5 into a vial of One Shot[®] Chemically Competent *E. coli* and mix gently. **Do not mix by pipetting up and down.**
 2. Incubate on ice for 5 minutes.
 3. Spread 50 μ l of cells on a prewarmed LB plate containing 50-100 μ g/ml ampicillin and incubate overnight at 37°C.
 4. An efficient TOPO[®] Cloning reaction should produce several hundred colonies. Pick ~10 white or light blue colonies for analysis (see Analyzing Positive Clones, page 12). Do not pick dark blue colonies.
-

continued on next page

Transforming One Shot[®] DH5 α [™] -T1^R, TOP10, and TOP10F['] Competent Cells, continued

One Shot[®] Electroporation Protocol

1. Add 2 μ l of the TOPO[®] Cloning reaction from Performing the TOPO[®] Cloning Reaction, Step 2, page 5 into a vial of One Shot[®] Electrocompetent *E. coli* and mix gently. **Do not mix by pipetting up and down.**
2. Carefully transfer solution to a 0.1 cm cuvette to avoid formation of bubbles.
3. Electroporate your samples using your own protocol and your electroporator.
Note: If you have problems with arcing, see below.
4. Immediately add 250 μ l of room temperature S.O.C. medium.
5. Transfer the solution to a 15 ml snap-cap tube (e.g. Falcon) and shake for at least 1 hour at 37°C to allow expression of the antibiotic resistance genes.
6. Spread 10-50 μ l from each transformation on a prewarmed selective plate and incubate overnight at 37°C. To ensure even spreading of small volumes, add 20 μ l of S.O.C. medium. We recommend that you plate two different volumes to ensure that at least one plate will have well-spaced colonies.
7. An efficient TOPO[®] Cloning reaction should produce several hundred colonies. Pick ~10 white or light blue colonies for analysis (see Analyzing Positive Clones, next page). Do not pick dark blue colonies.



Note

Addition of the Dilute Salt Solution in the TOPO[®] Cloning Reaction brings the final concentration of NaCl and MgCl₂ in the TOPO[®] Cloning reaction to 50 mM and 2.5 mM, respectively. To prevent arcing of your samples during electroporation, the volume of cells should be between 50 and 80 μ l (0.1 cm cuvettes) or 100 to 200 μ l (0.2 cm cuvettes).

If you experience arcing, try **one** of the following suggestions:

- Reduce the voltage normally used to charge your electroporator by 10%
- Reduce the pulse length by reducing the load resistance to 100 ohms
- Precipitate the TOPO[®] Cloning reaction and resuspend in water prior to electroporation

Analyzing Transformants

Analyzing Positive Clones

1. Take the 10 white or light blue colonies and culture them overnight in LB medium containing 50 µg/ml ampicillin or 50 µg/ml kanamycin.
Note: If you transformed One Shot[®] Mach1[™]-T1^R competent *E. coli*, you may inoculate overnight-grown colonies and culture them for 4 hours in **prewarmed** LB medium containing 50 µg/ml ampicillin or 50 µg/ml kanamycin before isolating plasmid. For optimal results, we recommend inoculating as much of a single colony as possible.
 2. Isolate plasmid DNA using your method of choice. If you need ultra-pure plasmid DNA for automated or manual sequencing, we recommend the S.N.A.P.[™] MiniPrep Kit (Catalog no. K1900-01) or the S.N.A.P.[™] MidiPrep Kit (Catalog no. K1910-01).
 3. Analyze the plasmids by restriction analysis to confirm the presence and correct orientation of the insert. Use a restriction enzyme or a combination of enzymes that cut once in the vector and once in the insert.
-

Sequencing

You may sequence your construct to confirm that your gene is cloned in the correct orientation. The M13 Forward (-20) and M13 Reverse primers are included to help you sequence your insert. Refer to the maps on page 15 (pCR[®]2.1-TOPO[®]) or page 16 (pCR[®]II-TOPO[®]) for sequence surrounding the TOPO TA Cloning[®] site. For the full sequence of either vector, refer to our Web site (www.invitrogen.com) or contact Technical Service (page 24).

Analyzing Transformants by PCR

You may wish to use PCR to directly analyze positive transformants. For PCR primers, use either the M13 Forward (-20) or the M13 Reverse primer and a primer that hybridizes within your insert. If you are using this technique for the first time, we recommend performing restriction analysis in parallel. Artifacts may be obtained because of mispriming or contaminating template. The protocol is provided below for your convenience. Other protocols are suitable.

Materials Needed

PCR SuperMix High Fidelity (Invitrogen, Catalog no. 10790-020)

Appropriate forward and reverse PCR primers (20 µM each)

Procedure

1. For each sample, aliquot 48 µl of PCR SuperMix High Fidelity into a 0.5 ml microcentrifuge tube. Add 1 µl each of the forward and reverse PCR primer.
 2. Pick 10 colonies and resuspend them individually in 50 µl of the PCR cocktail from Step 1, above. Don't forget to make a patch plate to preserve the colonies for further analysis.
 3. Incubate the reaction for 10 minutes at 94°C to lyse the cells and inactivate nucleases.
 4. Amplify for 20 to 30 cycles.
 5. For the final extension, incubate at 72°C for 10 minutes. Store at +4°C.
 6. Visualize by agarose gel electrophoresis.
-

continued on next page

Analyzing Transformants, continued



Important

If you have problems obtaining transformants or the correct insert, perform the control reactions described on pages 17-18. These reactions will help you troubleshoot your experiment.

Long-Term Storage

Once you have identified the correct clone, be sure to prepare a glycerol stock for long term storage. We recommend that you store a stock of plasmid DNA at -20°C .

1. Streak the original colony out on LB plates containing $50\ \mu\text{g/ml}$ ampicillin or $50\ \mu\text{g/ml}$ kanamycin.
 2. Isolate a single colony and inoculate into 1-2 ml of LB containing $50\ \mu\text{g/ml}$ ampicillin or kanamycin.
 3. Grow until culture reaches stationary phase.
 4. Mix 0.85 ml of culture with 0.15 ml of sterile glycerol and transfer to a cryovial.
 5. Store at -80°C .
-

Optimizing the TOPO[®] Cloning Reaction

Introduction

The information below will help you optimize the TOPO[®] Cloning reaction for your particular needs.

Faster Subcloning

The high efficiency of TOPO[®] Cloning technology allows you to streamline the cloning process. If you routinely clone PCR products and wish to speed up the process, consider the following:

- Incubate the TOPO[®] Cloning reaction for only 30 seconds instead of 5 minutes.

You may not obtain the highest number of colonies, but with the high efficiency of TOPO[®] Cloning, most of the transformants will contain your insert.

- After adding 2 μ l of the TOPO[®] Cloning reaction to chemically competent cells, incubate on ice for only 5 minutes.

Increasing the incubation time to 30 minutes does not significantly improve transformation efficiency.

More Transformants

If you are TOPO[®] Cloning large PCR products, toxic genes, or cloning a pool of PCR products, you may need more transformants to obtain the clones you want. To increase the number of colonies:

- Incubate the salt-supplemented TOPO[®] Cloning reaction for 20 to 30 minutes instead of 5 minutes.

Increasing the incubation time of the salt-supplemented TOPO[®] Cloning reaction allows more molecules to ligate, increasing the transformation efficiency. Addition of salt appears to prevent topoisomerase from rebinding and nicking the DNA after it has ligated the PCR product and dissociated from the DNA.

Cloning Dilute PCR Products



To clone dilute PCR products, you may:

- Increase the amount of the PCR product
 - Incubate the TOPO[®] Cloning reaction for 20 to 30 minutes
 - Concentrate the PCR product
-

Map of pCR[®]2.1-TOPO[®]

pCR[®]2.1-TOPO[®] Map

The map below shows the features of pCR[®]2.1-TOPO[®] and the sequence surrounding the TOPO[®] Cloning site. Restriction sites are labeled to indicate the actual cleavage site. The arrow indicates the start of transcription for T7 polymerase. **The complete sequence of pCR[®]2.1-TOPO[®] is available for downloading from our Web site (www.invitrogen.com) or by contacting Technical Service (page 24).**

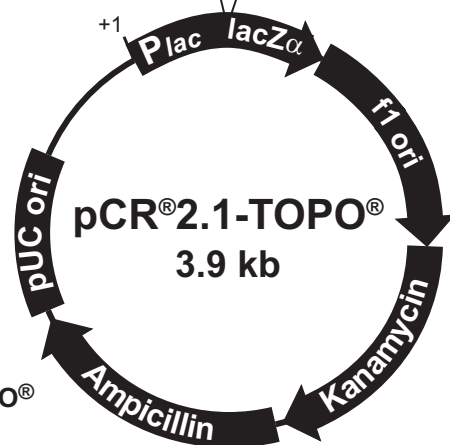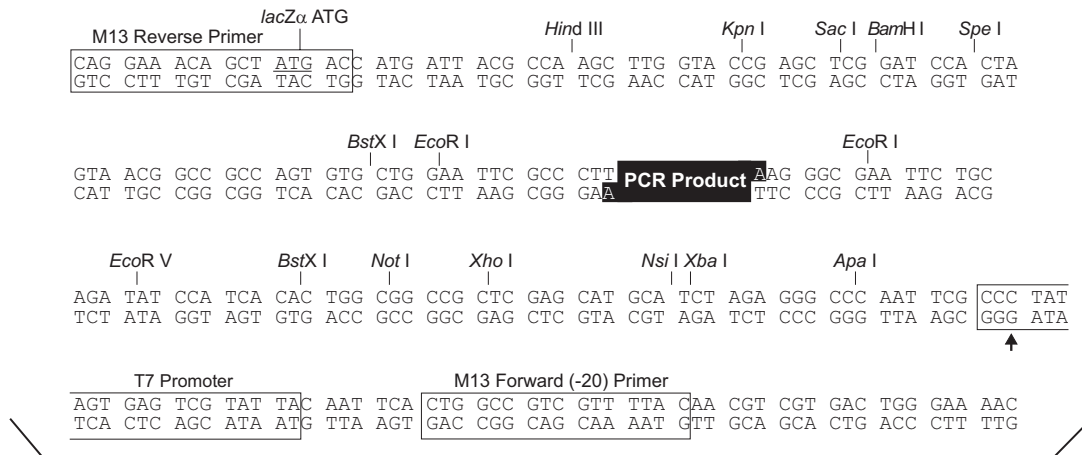


Comments for pCR[®]2.1-TOPO[®] 3931 nucleotides

- LacZα fragment: bases 1-547
- M13 reverse priming site: bases 205-221
- Multiple cloning site: bases 234-357
- T7 promoter/priming site: bases 364-383
- M13 Forward (-20) priming site: bases 391-406
- f1 origin: bases 548-985
- Kanamycin resistance ORF: bases 1319-2113
- Ampicillin resistance ORF: bases 2131-2991
- pUC origin: bases 3136-3809

Map of pCR[®]II-TOPO[®]

pCR[®]II-TOPO[®] Map The map below shows the features of pCR[®]II-TOPO[®] and the sequence surrounding the TOPO[®] Cloning site. Restriction sites are labeled to indicate the actual cleavage site. The arrows indicate the start of transcription for Sp6 and T7 polymerases. **The complete sequence of pCR[®]II-TOPO[®] is available for downloading from our Web site (www.invitrogen.com) or by contacting Technical Service (page 24).**

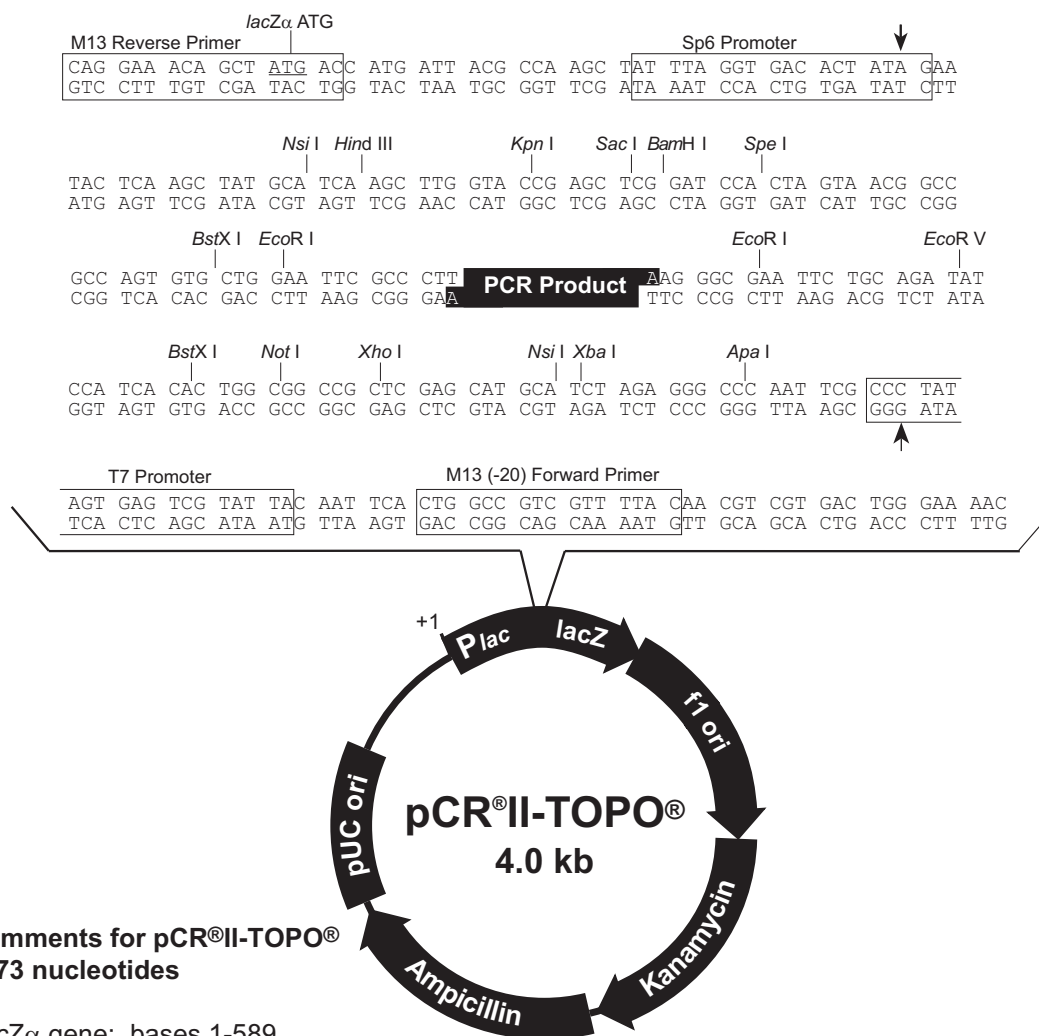

- LacZα* gene: bases 1-589
- M13 Reverse priming site: bases 205-221
- Sp6 promoter: bases 239-256
- Multiple Cloning Site: bases 269-383
- T7 promoter: bases 406-425
- M13 (-20) Forward priming site: bases 433-448
- f1 origin: bases 590-1027
- Kanamycin resistance ORF: bases 1361-2155
- Ampicillin resistance ORF: bases 2173-3033
- pUC origin: bases 3178-3851

Performing the Control Reactions

Introduction

We recommend performing the following control TOPO® Cloning reactions the first time you use the kit to help you evaluate your results. Performing the control reactions involves producing a control PCR product using the reagents included in the kit and using the PCR product directly in a TOPO® Cloning reaction.

Before Starting

For each transformation, prepare two LB plates containing 50 µg/ml kanamycin.

Note: Do not use plates containing ampicillin. The control template is a plasmid that encodes ampicillin resistance. This template is carried over into the TOPO® Cloning and transformation reactions. Transformants carrying this plasmid will also be ampicillin resistant and white, resulting in an apparent increase in TOPO® Cloning efficiency, but upon analysis, colonies do not contain the desired construct.

Producing Control PCR Product

- To produce the 750 bp control PCR product, set up the following 50 µl PCR:

Control DNA Template (100 ng)	1 µl
10X PCR Buffer	5 µl
dNTP Mix	0.5 µl
Control PCR Primers (0.1 µg/µl each)	1 µl
Sterile Water	41.5 µl
<u>Taq Polymerase (1 unit/µl)</u>	<u>1 µl</u>
Total Volume	50 µl
- Overlay with 70 µl (1 drop) of mineral oil.
- Amplify using the following cycling parameters:

Step	Time	Temperature	Cycles
Initial Denaturation	2 minute	94°C	1X
Denaturation	1 minute	94°C	25X
Annealing	1 minute	55°C	
Extension	1 minute	72°C	
Final Extension	7 minutes	72°C	1X

- Remove 10 µl from the reaction and analyze by agarose gel electrophoresis. A discrete 750 bp band should be visible. Proceed to the **Control TOPO® Cloning Reactions**, next page.
-

continued on next page

Performing the Control Reactions, continued

Control TOPO[®] Cloning Reactions

Using the control PCR product produced on the previous page and the TOPO[®] vector set up two 6 μ l TOPO[®] Cloning reactions as described below.

1. Set up control TOPO[®] Cloning reactions:

Reagent	"Vector Only"	"Vector + PCR Insert"
Control PCR Product	--	1 μ l
Sterile Water	4 μ l	3 μ l
Salt Solution or Dilute Salt Solution	1 μ l	1 μ l
TOPO [®] vector	1 μ l	1 μ l

2. Incubate at room temperature for **5 minutes** and place on ice.
3. Transform 2 μ l of each reaction into separate vials of One Shot[®] competent cells (pages 6-10).
4. Spread 10-50 μ l of each transformation mix onto LB plates containing 50 μ g/ml kanamycin and X-Gal (and IPTG, if using TOP10F' cells). Be sure to plate two different volumes to ensure that at least one plate has well-spaced colonies. For plating small volumes, add 20 μ l of S.O.C. medium to allow even spreading.
5. Incubate overnight at 37°C.

Analysis of Results

Hundreds of colonies from the vector + PCR insert reaction should be produced. Greater than ninety-five percent of these colonies will be white and 90% (or more) of these will contain the 750 bp insert when analyzed by EcoR I digestion and agarose gel electrophoresis.

Relatively few colonies will be produced in the vector-only reaction and most of these will be dark blue. You may observe a few white colonies. This results from removal of the 3' deoxythymidine overhangs creating a blunt-end vector. Ligation (re-joining) of the blunt ends will result in disruption of the *LacZ α* reading frame leading to the production of white colonies.

Transformation Control

pUC19 plasmid is included to check the transformation efficiency of the One Shot[®] competent cells. Transform with 10 pg per 50 μ l of cells using the protocols on pages 6-10.

Use LB plates containing 100 μ g/ml ampicillin. Just before plating the transformation mix for electrocompetent cells, dilute 10 μ l of the mix with 90 μ l of S.O.C. medium.

Type of Cells	Volume to Plate	Transformation Efficiency
Chemically Competent	10 μ l + 20 μ l S.O.C.	$\sim 1 \times 10^9$ cfu/ μ g DNA
Electrocompetent	20 μ l (1:10 dilution)	$> 1 \times 10^9$ cfu/ μ g DNA

continued on next page

Performing the Control Reactions, continued

Factors Affecting Cloning Efficiency

Note that lower cloning efficiencies will result from the following variables. Most of these are easily correctable, but if you are cloning large inserts, you may not obtain the expected 95% (or more) cloning efficiency.

Variable	Solution
pH>9	Check the pH of the PCR amplification reaction and adjust with 1 M Tris-HCl, pH 8.
Incomplete extension during PCR	Be sure to include a final extension step of 7 to 30 minutes during PCR. Longer PCR products will need a longer extension time.
Cloning large inserts (>1 kb)	Try one or all of the following: Increase amount of insert. Incubate the TOPO [®] Cloning reaction longer. Gel-purify the insert (see page 20).
Excess (or overly dilute) PCR product	Reduce (or concentrate) the amount of PCR product.
Cloning blunt-ended fragments	Add 3' A-overhangs to your blunt PCR product by incubating with <i>Taq</i> polymerase (page 22).
	Use the Zero Blunt [®] PCR Cloning Kit to clone blunt PCR products (Catalog no. K2700-20).
PCR cloning artifacts ("false positives")	TOPO [®] Cloning is very efficient for small fragments (< 100 bp) present in certain PCR reactions. Gel-purify your PCR product (page 20).
PCR product does not contain sufficient 3' A-overhangs even though you used <i>Taq</i> polymerase	Increase the final extension time to ensure all 3' ends are adenylated. <i>Taq</i> polymerase is less efficient at adding a nontemplate 3' A next to another A. <i>Taq</i> is most efficient at adding a nontemplate 3' A next to a C. You may have to redesign your primers so that they contain a 5' G instead of a 5' T (Brownstein <i>et al.</i> , 1996).

Appendix

Purifying PCR Products

Introduction

Smearing, multiple banding, primer-dimer artifacts, or large PCR products (>1 kb) may necessitate gel purification. If you intend to purify your PCR product, be extremely careful to remove all sources of nuclease contamination. There are many protocols to isolate DNA fragments or remove oligonucleotides. Refer to Current Protocols in Molecular Biology, Unit 2.6 (Ausubel et al., 1994) for the most common protocols. Two simple protocols are provided below for your convenience.

Using the S.N.A.P.[™] MiniPrep Kit

The S.N.A.P.[™] MiniPrep Kit (Catalog no. K1900-01) allows you to rapidly purify PCR products from regular agarose gels. You will need to prepare 6 M sodium iodide, 10 mM sodium sulfite solution in sterile water before starting. Sodium sulfite prevents oxidation of NaI.

1. Electrophorese amplification reaction on a 1 to 5% regular TAE agarose gel.
Note: Do not use TBE. Borate interferes with the NaI step (Step 2).
 2. Cut out the gel slice containing the PCR product and melt it at 65°C in 2 volumes of 6 M NaI.
 3. Add 1.5 volumes Binding Buffer (provided in the S.N.A.P.[™] MiniPrep Kit).
 4. Load solution (no more than 1 ml at a time) from Step 3 onto a S.N.A.P.[™] column. Centrifuge 1 minute at full speed in a microcentrifuge and discard the supernatant.
 5. If you have solution remaining from Step 3, repeat Step 4.
 6. Add 900 µl of the Final Wash Buffer (provided in the S.N.A.P.[™] MiniPrep Kit).
 7. Centrifuge 1 minute at full speed in a microcentrifuge and discard the supernatant. Repeat.
 8. Elute the purified PCR product in 40 µl of TE or sterile water. Use 4 µl for the TOPO[®] Cloning reaction and proceed as described on page 5.
-

Quick S.N.A.P.[™] Method

An easier method is to simply cut out the gel slice containing your PCR product, place it on top of the S.N.A.P.[™] column bed, and centrifuge at full speed for 10 seconds. Use 1-2 µl of the flow-through in the TOPO[®] Cloning reaction (page 5). Be sure to make the gel slice as small as possible for best results.

continued on next page

Purifying PCR Products, continued

Low-Melt Agarose Method

Note that gel purification will result in a dilution of your PCR product. Use only chemically competent cells for transformation.

1. Electrophorese all of your PCR reaction on a low-melt TAE agarose gel (0.8 to 1.2%).
 2. Visualize the band of interest and excise the band.
 3. Place the gel slice in a microcentrifuge tube and incubate the tube at 65°C until the gel slice melts.
 4. Place the tube at 37°C to keep the agarose melted.
 5. Use 4 μ l of the melted agarose containing your PCR product in the TOPO[®] Cloning reaction (page 5).
 6. Incubate the TOPO[®] Cloning reaction **at 37°C for 5 to 10 minutes**. This is to keep the agarose melted.
 7. Transform 2 to 4 μ l directly into competent One Shot[®] cells using one of the methods described on pages 10-11.
-



Note

Note that the cloning efficiency may decrease with purification of the PCR product. You may wish to optimize your PCR to produce a single band.

Addition of 3' A-Overhangs Post-Amplification

Introduction

Direct cloning of DNA amplified by proofreading polymerases into TOPO TA Cloning[®] vectors is often difficult because proofreading polymerases remove the 3' A-overhangs necessary for TA Cloning[®]. Invitrogen has developed a simple method to clone these blunt-ended fragments.

Before Starting

You will need the following items:

- *Taq* polymerase
 - A heat block equilibrated to 72°C
 - Phenol-chloroform (optional)
 - 3 M sodium acetate (optional)
 - 100% ethanol (optional)
 - 80% ethanol (optional)
 - TE buffer (optional)
-

Procedure

This is just one method for adding 3' adenines. Other protocols may be suitable.

1. After amplification with a proofreading polymerase, place vials on ice and add 0.7-1 unit of *Taq* polymerase per tube. Mix well. It is not necessary to change the buffer. A sufficient number of PCR products will retain the 3' A-overhangs.
2. Incubate at 72°C for 8-10 minutes (do not cycle).
3. Place on ice and use immediately in the TOPO[®] Cloning reaction.

Note: If you plan to store your sample overnight before proceeding with TOPO[®] Cloning, extract your sample with an equal volume of phenol-chloroform to remove the polymerases. Ethanol-precipitate the DNA and resuspend in TE buffer using the starting volume of the PCR.



Note

You may also gel-purify your PCR product after amplification with a proofreading polymerase. After purification, add *Taq* polymerase buffer, dATP, and 0.5 unit of *Taq* polymerase. Incubate the reaction for 10-15 minutes at 72°C and use in the TOPO[®] Cloning reaction.

Recipes

LB (Luria-Bertani) Medium and Plates

Composition:

1.0% Tryptone
0.5% Yeast Extract
1.0% NaCl
pH 7.0

1. For 1 liter, dissolve 10 g tryptone, 5 g yeast extract, and 10 g NaCl in 950 ml deionized water.
2. Adjust the pH of the solution to 7.0 with NaOH and bring the volume up to 1 liter.
3. Autoclave on liquid cycle for 20 minutes at 15 psi. Allow solution to cool to 55°C and add antibiotic if needed (50 µg/ml of either ampicillin or kanamycin).
4. Store at room temperature or at +4°C.

LB agar plates

1. Prepare LB medium as above, but add 15 g/L agar before autoclaving.
 2. Autoclave on liquid cycle for 20 minutes at 15 psi.
 3. After autoclaving, cool to ~55°C, add antibiotic (50 µg/ml of either ampicillin or kanamycin), and pour into 10 cm plates.
 4. Let harden, then invert and store at +4°C in the dark.
-

Technical Service

World Wide Web



Visit the Invitrogen Web Resource using your World Wide Web browser. At the site, you can:

- Get the scoop on our hot new products and special product offers
- View and download vector maps and sequences
- Download manuals in Adobe® Acrobat® (PDF) format
- Explore our catalog with full color graphics
- Obtain citations for Invitrogen products
- Request catalog and product literature

Once connected to the Internet, launch your Web browser (Internet Explorer 5.0 or newer or Netscape 4.0 or newer), then enter the following location (or URL):

<http://www.invitrogen.com>

...and the program will connect directly. Click on underlined text or outlined graphics to explore. Don't forget to put a bookmark at our site for easy reference!

Contact Us

For more information or technical assistance, call, write, fax, or email. Additional international offices are listed on our Web page (www.invitrogen.com).

Corporate Headquarters:

Invitrogen Corporation
1600 Faraday Avenue
Carlsbad, CA 92008 USA
Tel: 1 760 603 7200
Tel (Toll Free): 1 800 955 6288
Fax: 1 760 602 6500
E-mail:
tech_service@invitrogen.com

Japanese Headquarters:

Invitrogen Japan K.K.
Nihonbashi Hama-Cho Park Bldg. 4F
2-35-4, Hama-Cho, Nihonbashi
Tel: 81 3 3663 7972
Fax: 81 3 3663 8242
E-mail: jpinfo@invitrogen.com

European Headquarters:

Invitrogen Ltd
Inchinnan Business Park
3 Fountain Drive
Paisley PA4 9RF, UK
Tel: +44 (0) 141 814 6100
Tech Fax: +44 (0) 141 814 6117
E-mail: eurotech@invitrogen.com

MSDS Requests

To request an MSDS, visit our Web site at www.invitrogen.com. On the home page, go to 'Technical Resources', select 'MSDS', and follow instructions on the page.

continued on next page

Technical Service, continued

Limited Warranty

Invitrogen is committed to providing our customers with high-quality goods and services. Our goal is to ensure that every customer is 100% satisfied with our products and our service. If you should have any questions or concerns about an Invitrogen product or service, please contact our Technical Service Representatives.

Invitrogen warrants that all of its products will perform according to the specifications stated on the certificate of analysis. The company will replace, free of charge, any product that does not meet those specifications. This warranty limits Invitrogen Corporation's liability only to the cost of the product. No warranty is granted for products beyond their listed expiration date. No warranty is applicable unless all product components are stored in accordance with instructions. Invitrogen reserves the right to select the method(s) used to analyze a product unless Invitrogen agrees to a specified method in writing prior to acceptance of the order.

Invitrogen makes every effort to ensure the accuracy of its publications, but realizes that the occasional typographical or other error is inevitable. Therefore Invitrogen makes no warranty of any kind regarding the contents of any publications or documentation. If you discover an error in any of our publications, please report it to our Technical Service Representatives.

Invitrogen assumes no responsibility or liability for any special, incidental, indirect or consequential loss or damage whatsoever. The above limited warranty is sole and exclusive. No other warranty is made, whether expressed or implied, including any warranty of merchantability or fitness for a particular purpose.

Product Qualification

Restriction Digest Supercoiled pCR[®]2.1-TOPO[®] and pCR[®]II-TOPO[®] are qualified by restriction digest. The table below lists the restriction enzymes and the expected fragments.

Restriction Enzyme	pCR [®] 2.1-TOPO [®]	pCR [®] II-TOPO [®]
<i>Hind</i> III (linearizes)	3890 bp	3932 bp
<i>Xba</i> I (linearizes)	3890 bp	3932 bp
<i>Nsi</i> I	3890 bp	96, 3836 bp
<i>Pst</i> I	1167, 2723 bp	1167, 2765 bp
<i>EcoR</i> I and <i>Afl</i> III	408, 693, 2789bp	450, 693, 2789 bp

TOPO[®] Cloning Efficiency

Once the vectors have been adapted with topoisomerase I, they are lot-qualified using the control reagents included in the kit. Under conditions described on pages 17-18, a 750 bp control PCR product was TOPO[®]-Cloned into each vector and subsequently transformed into the One Shot[®] competent *E. coli* included with the kit.

Each lot of vector should yield greater than 95% cloning efficiency.

Primers

Both primers have been lot-qualified by DNA sequencing experiments using the dideoxy chain termination technique.

One Shot[®] Competent *E. coli*

All competent cells are qualified as follows:

- Cells are tested for transformation efficiency using the control plasmid included in the kit. Transformed cultures are plated on LB plates containing 100 µg/ml ampicillin and the transformation efficiency is calculated. Test transformations are performed in duplicate. Transformation efficiency should be $\sim 1 \times 10^9$ cfu/µg DNA for chemically competent cells and $>1 \times 10^9$ for electrocompetent cells.
 - To verify the absence of phage contamination, 0.5-1 ml of competent cells are added to LB top agar and poured onto LB plates. After overnight incubation, no plaques should be detected.
 - Untransformed cells are plated on LB plates 100 µg/ml ampicillin, 25 µg/ml streptomycin, 50 µg/ml kanamycin, or 15 µg/ml chloramphenicol to verify the absence of antibiotic-resistant contamination.
-

Purchaser Notification

Information for European Customers

The Mach1™-T1^R *E. coli* strain is genetically modified to carry the *lacZ*ΔM15 *hsdR lacX74 recA endA tonA* genotype. As a condition of sale, this product must be in accordance with all applicable local legislation and guidelines including EC Directive 90/219/EEC on the contained use of genetically modified organisms.

Limited Use Label License No. 118: TOPO® Cloning Products

The TOPO® Cloning Technology products and their use are the subject of one or more of U.S. Patent Nos. 5,766,891, 6,548,277 and/or other pending U.S. and foreign patent applications licensed to Invitrogen Corporation.

The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. The buyer may transfer information or materials made through the use of this product to a scientific collaborator, provided that such transfer is not for any Commercial Purpose, and that such collaborator agrees in writing (a) not to transfer such materials to any third party, and (b) to use such transferred materials and/or information solely for research and not for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. Invitrogen Corporation will not assert a claim against the buyer of infringement of the above patents based upon the manufacture, use or sale of a therapeutic, clinical diagnostic, vaccine or prophylactic product developed in research by the buyer in which this product or its components was employed, provided that neither this product nor any of its components was used in the manufacture of such product. If the purchaser is not willing to accept the limitations of this limited use statement, Invitrogen is willing to accept return of the product with a full refund. For information on purchasing a license to this product for purposes other than research, contact Licensing Department, Invitrogen Corporation, 1600 Faraday Avenue, Carlsbad, California 92008. Phone (760) 603-7200. Fax (760) 602-6500.

continued on next page

Purchaser Notification, continued

**Limited Use Label
License No. 120:
TA[®] Cloning
Products**

The TA[®] Cloning Technology products and their use are the subject of one or more of U.S. Patent Nos. 5,487,993, and 5,827,657, and/or pending foreign patent applications owned by Invitrogen Corporation.

The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. The buyer may transfer information or materials made through the use of this product to a scientific collaborator, provided that such transfer is not for any Commercial Purpose, and that such collaborator agrees in writing (a) not to transfer such materials to any third party, and (b) to use such transferred materials and/or information solely for research and not for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. Invitrogen Corporation will not assert a claim against the buyer of infringement of the above patents based upon the manufacture, use or sale of a therapeutic, clinical diagnostic, vaccine or prophylactic product developed in research by the buyer in which this product or its components was employed, provided that neither this product nor any of its components was used in the manufacture of such product. If the purchaser is not willing to accept the limitations of this limited use statement, Invitrogen is willing to accept return of the product with a full refund. For information on purchasing a license to this product for purposes other than research, contact Licensing Department, Invitrogen Corporation, 1600 Faraday Avenue, Carlsbad, California 92008. Phone (760) 603-7200. Fax (760) 602-6500.

continued on next page

Purchaser Notification, continued

**Limited Use Label
License No. 144:
Mach1™
Competent Cells**

The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for commercial purposes. The buyer may transfer information or materials made through the use of this product to a scientific collaborator, provided that such transfer is not for any commercial purpose, and that such collaborator agrees in writing (a) to not transfer such materials to any third party, and (b) to use such transferred materials and/or information solely for research and not for commercial purposes. Commercial purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. Invitrogen Corporation will not assert a claim against the buyer of infringement of patents owned by Invitrogen and claimed in this product based upon the manufacture, use or sale of a therapeutic, clinical diagnostic, vaccine or prophylactic product developed in research by the buyer in which this product or its components was employed, provided that neither this product nor any of its components was used in the manufacture of such product. If the purchaser is not willing to accept the limitations of this limited use statement, Invitrogen is willing to accept return of the product with a full refund. For information on purchasing a license to this product for purposes other than research, contact Licensing Department, Invitrogen Corporation, 1600 Faraday Avenue, Carlsbad, California 92008. Phone (760) 603-7200. Fax (760) 602-6500.

References

- Ausubel, F. M., Brent, R., Kingston, R. E., Moore, D. D., Seidman, J. G., Smith, J. A., and Struhl, K. (1994). *Current Protocols in Molecular Biology* (New York: Greene Publishing Associates and Wiley-Interscience).
- Brownstein, M. J., Carpten, J. D., and Smith, J. R. (1996). Modulation of Non-Templated Nucleotide Addition by Taq DNA Polymerase: Primer Modifications that Facilitate Genotyping. *BioTechniques* 20, 1004-1010.
- Innis, M. A., Gelfand, D. H., Sninsky, J. J., and White, T. S. (1990) *PCR Protocols: A Guide to Methods and Applications*. Academic Press, San Diego, CA.
- Sambrook, J., Fritsch, E. F., and Maniatis, T. (1989). *Molecular Cloning: A Laboratory Manual*, Second Edition (Plainview, New York: Cold Spring Harbor Laboratory Press).
- Shuman, S. (1991). Recombination Mediated by Vaccinia Virus DNA Topoisomerase I in *Escherichia coli* is Sequence Specific. *Proc. Natl. Acad. Sci. USA* 88, 10104-10108.
- Shuman, S. (1994). Novel Approach to Molecular Cloning and Polynucleotide Synthesis Using Vaccinia DNA Topoisomerase. *J. Biol. Chem.* 269, 32678-32684.

©1999-2004 Invitrogen Corporation. All rights reserved.

For research use only. Not intended for any animal or human therapeutic or diagnostic use.

pCR[®], TOPO[®], and TOPO TA Cloning[®] are registered trademarks of Invitrogen Corporation.